

## **IMPORTANT SAFEGUARDS**

All information in this section is provided to instruct you in the correct and safe operation of this device so that you can avoid injuries to yourself and others as well as damage to property.

### **READ ALL INSTRUCTIONS BEFORE USING**

#### **DANGERS**


- Toddlers, infants, people who are unable to operate should not allowed use this appliance. Please keep it out of reach of children.
- Never alter the appliance. Also, do not try to disassemble or repair it.
- Do not place or store the product where it can fall or be pulled into tub or sink.
- Do not reach for the product that has fallen into water. Unplug immediately.
- Do not use the product while bathing.

#### **WARNINGS**

- Oral irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- Children shall not play with the appliance.
- The appliance is only to be used with the power supply unit provided with the appliance.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery;
- The battery is to be disposed of safely.
- Regarding how to remove the battery, please to see “Disposal” in page 4.
- Do not fill the tank with water warmer than 40 °C
- Never use the adaptor if it has a damaged cord or plug.
- Keep the cord away from heated surfaces.
- Please do not use it if you lost the nozzle.
- Do not submerge the device in water.
- Do not use for any purpose other than oral cleaning.
- Keep nozzle away from small children. The nozzle is small enough to be swallowed and can become lodged in the throat.
- Certain types of mouthwash may damage the equipment by causing the housing or water reservoir to crack. Therefore, do not use mouthwash in the water reservoir.
- People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- Do not press the nozzle too firmly against your teeth or gums.
- Doing so may injure your gums.

### **SAVE THESE INSTRUCTIONS!**

|   |                       |
|---|-----------------------|
| A | Nozzle                |
| B | Nozzle Release Button |
| C | ON/OFF Button         |
| D | Mode Button           |
| E | Mode Indicator        |
| F | Water Tank            |
| G | Charging Base Cover   |
| H | Charging Base         |
| I | Water Tank Cup        |

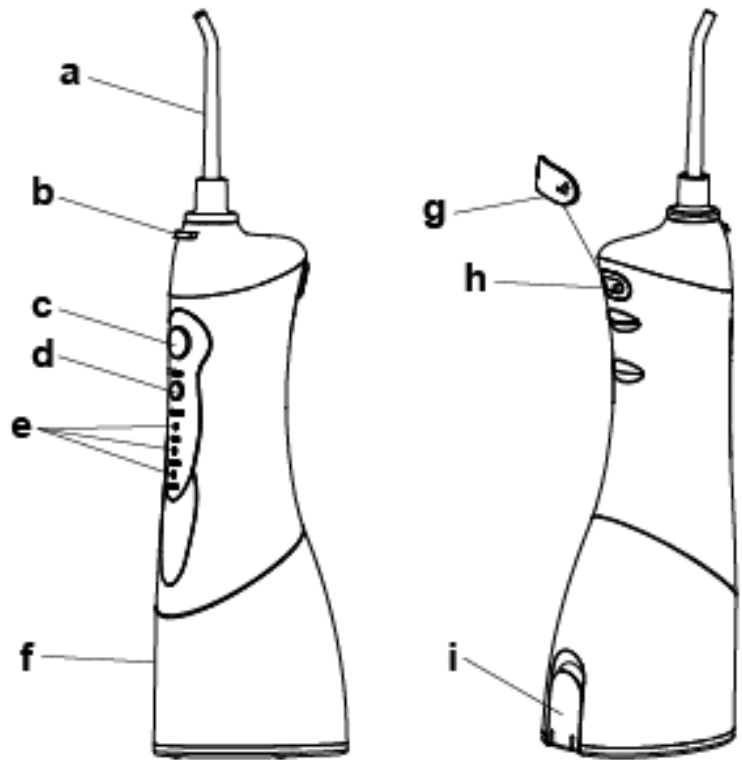


**Note:**

There may be some liquid remaining inside the main appliance before use (tank, water suction tube, etc), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

**How to use**

**Charging the device**

1. Take the charger base cover down and insert the appliance plug securely into the appliance socket (make sure no liquid remaining inside the socket before insert the appliance plug).
2. Insert the charging adaptor securely into the power outlet.
3. The “ch” indicator blinks when the device is charging, and it will be steadily on after device is full charged.
4. Recover the cover onto charger base again after charging finished

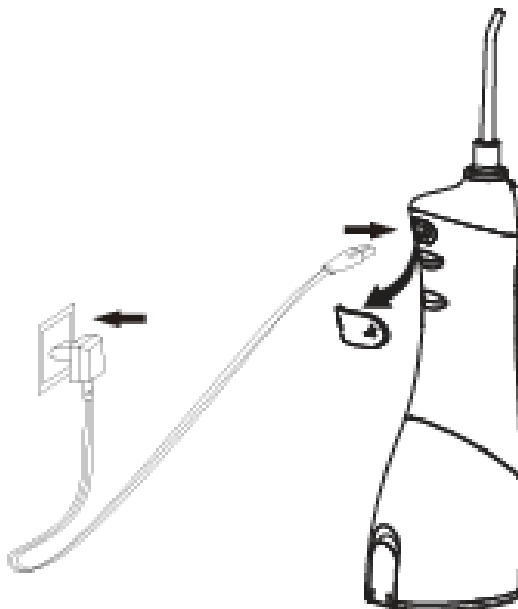


**Note:**


If the “CH” indicator is blinking during or after an irrigating, that means the battery need to be recharged. An full recharging will take about 8 hours. (Fig.1)

## Using the Water Flosser:


Attach nozzle onto the unit.

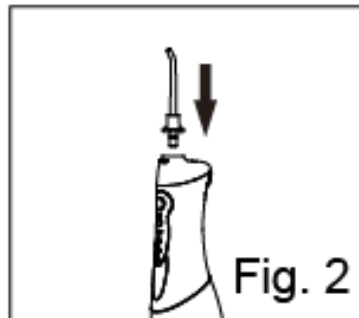


Open the water tank cap to fill with water, which temperature should not exceed 40°C, then close the cap.

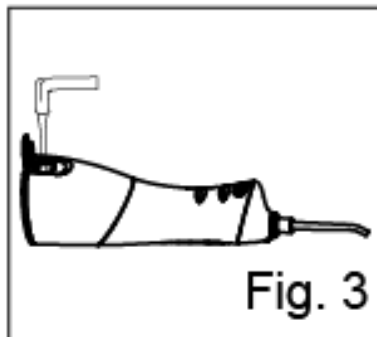

You can also slide the water tank off of the unit before refilling.

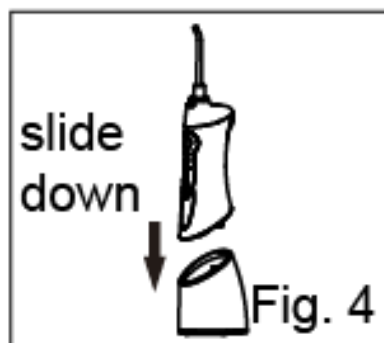

Press the “Mode” key to select between the desired flossing mode.

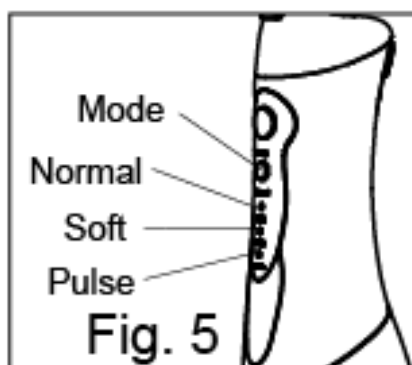

### Flossing Modes:

**Normal:** for normal cleaning.

**Soft:** for gentle cleaning of sensitive gums.


**Pulse:** for cleaning and massaging the gums.

## Using the Water Flosser:

Hold the unit in an upright position with nozzle tip point toward your teeth


Keep your mouth opened slightly so that the water can flow out.

Press “On/Off” key to start irrigating. (Fig. 6)

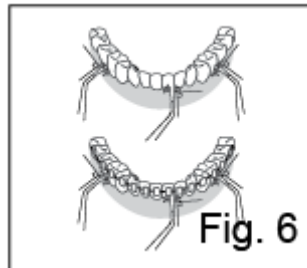


- 1). Directing the water jet flow at a 90-degree angle to the teeth and gums
- 2). Directing the water jet slowly along the teeth
- 3). Directing the water jet along the gum line with the nozzle tip at a 90-degree angle to the gum tissue
- 4). Cleaning between the teeth and around any braces, crowns or bridgework thoroughly.

Press “On/Off” key and then press the nozzle release button to replace other nozzle when you need.(Fig.7)

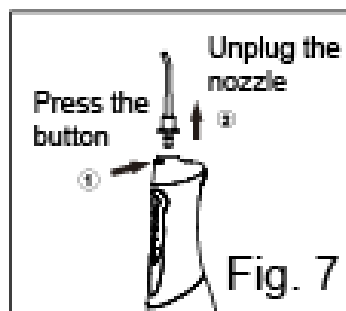

Once irrigating complete, press “On/Off” button again to turn off the device.

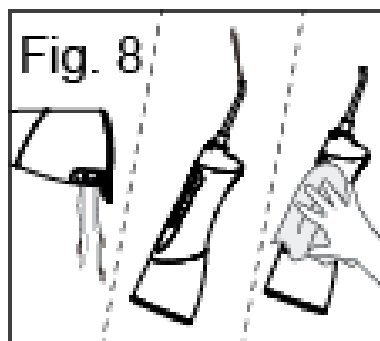

1. Open the water tank cap and dispose of the water in the tank.
2. Press “On/Off” key to turn on the unit, and drain the water inside the appliance.
3. Press “On/Off” key again to turn off the unit, and close the water tank cap.
4. Wipe the device with dry cloth.

**Note:**

To prevent bacteria propagating, make sure there is no remaining water inside the unit after use. If you're not going to use this appliance for long time, be sure wipe and dry the appliance immediately after use.

**Maintenance**

Use only water or neutral cleaning agents to clean the device. Never use caustic, abrasive/harsh agents/chemicals (e.g., vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 50°C. Keep the device away from high temperature areas and avoid direct sunlight.

**Main unit**

- 1.Remove the dirt by wiping using a cloth with liquid soap and wipe dry using a dry cloth.
- 2.Do not submerge the main unit in water to clean.
- 3.Dry after washing with water.

**Water tank cap**

Wipe off any grime and dirt with a cloth.

**Nozzle**

1. Wash with water and wipe with soft cloth.
2. Replace the nozzle every 6 months.

**Water suction tube**

1. Wash with water and wipe with soft cloth.
2. Do not bend, pull or twist the water suction tube.

**Tank**

1. Wash with water and wipe with soft cloth.
2. Dry the tank when not using for a long period.

**Troubleshooting**

| Problem   | Possible Cause  | Solution                                   |
|---|---|--|
| It does not work                                    | You've just purchased the product or it has been left unused for more than 3 months | Charge the unit for at least 8 hours.      |
| It only works for a few minutes even after charging | The battery has reached the end of its life cycle.                                  | Contact service center for support.        |
|   | Charging time was insufficient.   | Charge the unit for at least 8 hours.      |
| Water pressure is very low                          | There is dirt stuck in the nozzle   | Replace the nozzle                         |
|   | The nozzle is deformed  |  |
|   | There is dirt stuck in the filter   | Clean the filter                           |
| Water does not come out                             | The tank is empty   | Fill up tank                               |
|   | The appliance is tilted too much when using   | Use with the appliance in a vertical state |